




Yes, you can.™

INVACARE PORTABLE RAMPS



INVACARE® Portable Ramps

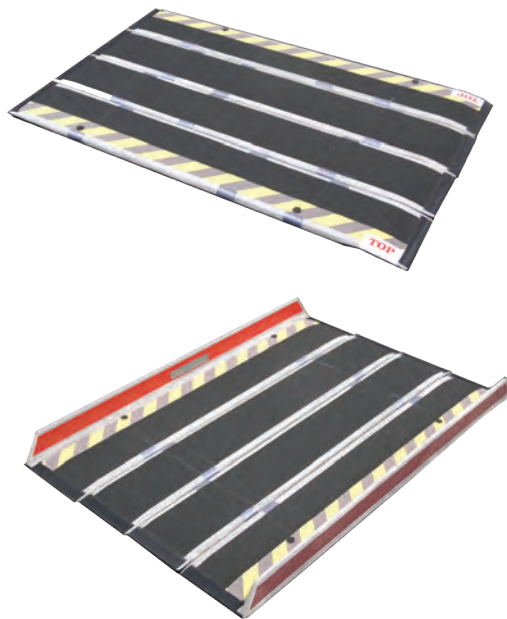


Invacare Portable Ramps are manufactured with a unique combination of uni-directional glass and carbon fibre materials to create an extremely strong and lightweight ramp suitable for a wide range of uses.


Tested to AS/NZ Standards, up to 300kg, Invacare Portable Ramps can be deployed in seconds, requiring no assembly and is immediately ready for use!*



Carry bag for Personal and Multipurpose models

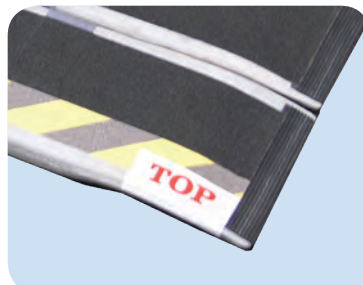

Marked 'TOP' to show head of ramp, grip tape, pads at both ends

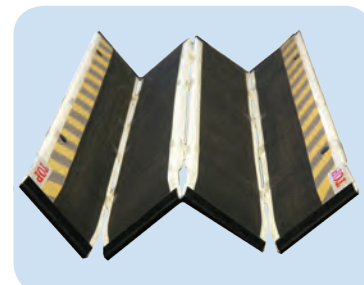

Multi-paneled, hinge system folds with a concertina action, allowing it to be folded to a single compact unit

Features

Invacare Portable Ramps can be quickly folded into a single compact unit.

Made from a unique combination of uni-directional glass and carbon fibre materials.

Multi-paneled design held together with a patented hinge system connected by webbing

Clearly marked and grip tape covered, Invacare Portable Ramps are also fitted with pads at the top and base.

Invacare Portable Ramps are tested to AS/NZ Standards, up to 300kg!*

Invacare Portable Ramps require no assembly.

Benefits

PORTABLE - Allowing for easy storage or transport in any vehicle or home environment.

LIGHTWEIGHT - Invacare Portable Ramps are very strong, yet light enough to be portable.

COMPACT - Enables Invacare Portable Ramps to fold with a concertina action, allowing for a much smaller width when folded.

SAFE - Invacare Portable Ramps provide a non-slip surface, for all weather conditions.

STRONG - For security and peace of mind.

SIMPLE-TO-USE - Allowing the ramp to be deployed in seconds, and immediately ready for use!



INVACARE®

Portable Ramps

SIMPLE-TO-USE - Invacare Portable Ramps can be deployed in seconds, requiring no assembly and is immediately ready for use!



Invacare Portable Ramps are available in 4 models, in a wide range of sizes;

PERSONAL - Portable, lightweight ramp that allows easier access in and around the home or workplace. Ideal for climbing kerbs or similar obstacles. Includes carry bag as standard.

MULTIPURPOSE - Longer and wider, designed to cope with a wider range of obstacles making it more suitable for commercial and travel situations. Invaluable when venturing outdoors. Includes carry bag as standard.

SENIOR - Flexible range of ramps to cope with obstacles from 22.5-87cm in height. 10 sizes available.

EDGE BARRIER LIMITER (EBL) - Added security of edge barriers, particularly suited to providing instant access to buses, trains and buildings where a permanent ramp is impractical. 9 sizes available.

Model/Length (cm)	Max Height Capacity	Ramp Width	Ramp Weight
Personal 070	17.5cm	74.5cm	4.7kg
Multipurpose 087	21.75cm	83.5cm	5.6kg
Multipurpose 100	25cm	83cm	6kg
Senior 090	22.5cm	75cm	4.8kg
Senior 120	30cm	75cm	8.3kg
Senior 135	33.75cm	75cm	8.5kg
Senior 150	37.5cm	75cm	9.6kg
Senior 165	41.2cm	75cm	10.8kg
Senior 185	46.25cm	75cm	11.2kg
Senior 200	50cm	75cm	11.7kg
Senior 250	62.5cm	75cm	14.8kg
Senior 300	75cm	75cm	22.1kg
Senior 350	87.5cm	75cm	24.7kg
EBL 090	22.5cm	75cm	6kg
EBL 120	30cm	82cm	9.7kg
EBL 135	33.75cm	82cm	11.1kg
EBL 150	37.5cm	82cm	11.8kg
EBL 165	41.2cm	82cm	13kg
EBL 200	50cm	82cm	13.7kg
EBL 250	62.5cm	82cm	23kg
EBL 300	75cm	82cm	27.3kg
EBL 350	87.5cm	82cm	31.5kg

**Nominal weights may vary due to local specification. Please confirm with your Invacare Portable Ramp Dealer.*

INVACARE®

COMPANY PROFILE

INVACARE CORPORATION traces its roots back to 1885 when the Worthington Company began producing a line of vehicles designed for people with physical disabilities. Through the years, the complexion of the company has changed several times and the product range has expanded dramatically.

Today, Invacare is the global leader in the \$6 billion market for home medical products, manufacturing and distributing the world's broadest product offering to approximately 15,000 independent, home medical equipment providers.

Since establishment in Australia, Invacare has grown both organically and through the acquisition of several trusted and well known Australian brands including Roller Chair, Pro-Med, Decpac, Australian Healthcare Equipment (Murphy Furniture and Hendicare Hospital Equipment), Morris Surgical, Home Health Equipment and Healthcare Equipment WA.

The acquisition of these Australian organisations has created a strong local business that has a broad product portfolio to support institutional and home healthcare.

Invacare Australia's current product range includes electric beds, furniture, static and active pressure care mattresses, patient lifters and slings, oxygen concentrators, manual and power wheelchairs, lightweight portable ramps, walking, toilet and bathroom aids.



Yes, you can.™

Invacare New Zealand
4 Westfield Place (PO Box 62-124) Mt Wellington Auckland, New Zealand
Freephone: 0508 468 222 Freefax: 0508 807 788
Email: sales@invacare.co.nz Web: www.invacare.co.nz

Invacare Australia Pty Ltd
1 Lenton Place (PO Box 5002) North Rocks NSW 2151, Australia
Freephone: 1800 460 460 Fax: 02 8839 5353
Email: sales@invacare.com.au Web: www.invacare.com.au

Invacare Asia Ltd
1 Lenton Place (PO Box 5002) North Rocks NSW 2151, Australia
Phone: +61 2 8839 5330 Fax: +61 2 8839 5331
Email: asiasales@invacare.com Web: www.invacare.com