

Hyper Foam 2


Decubitus therapy system

autoclavable according to
RKI*¹ guidelines

applicable for
pain therapy

up to stage 4
EPUAP

high risk

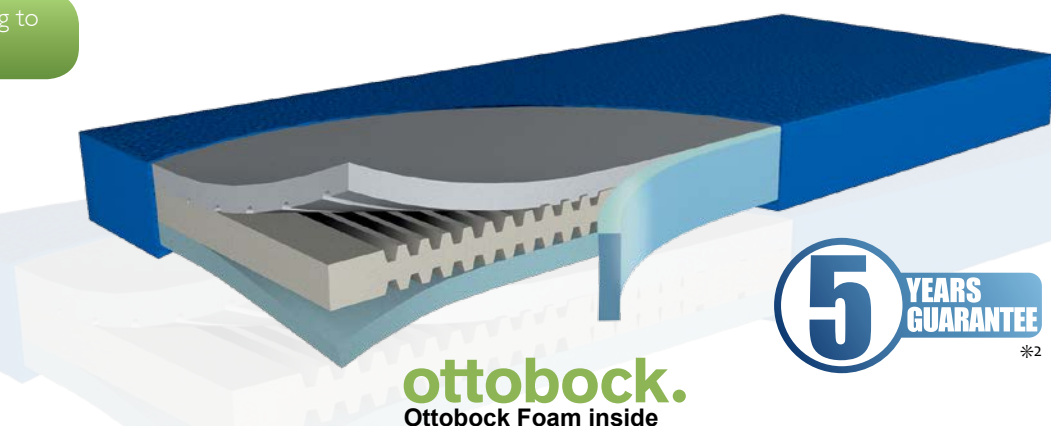

5 YEARS
GUARANTEE
*²

- hyperelastic HR-foam core with longitudinal and diagonal intersections
- three-layered
- exceptional ventilation system for an ideal microclimate
- optimum pressure relief and minimum shear forces (below 18mmHg pressure)
- PU-Tex cover removable and washable up to 95°C (autoclavable + disinfectant-proof)
- mattress autoclavable: chemical and thermal treatment according to the directives of the Robert-Koch-Institute possible
- permeable to water vapour
- flame retardant according to DIN EN 597 (part I + 2)/CRIB 5
- suitable for patients with a body weight from 20 kg to 120 kg
- applicable for pain therapy

For sizes have a look in the current price list.



*¹ Robert-Koch-Institute (RKI) is a part of the Federal Ministry of Health in Germany. It is responsible for disease prevention and control.

*² Guarantee applies to functional properties of the mattress core!